

Data Sheet – SFM3100 Low Pressure Drop Analog Flow Meter

- Analog output
- Compact design
- Flow range: -24 ... 240 slm
- Direct gas temperature measurement
- Fully calibrated & temperature compensated
- Very fast response time



Product Summary

The SFM3100 sensor is Sensirion's flow meter designed for the mass flow measurement in ventilators. It measures the flow rate of **air**, **oxygen and non-aggressive gases** with superb accuracy. A special design of the flow channel results in the very low pressure drop through the flow body of the sensor making it extremely suitable for medical ventilation and respiratory applications.

The SFM3100 operates from a 5 V supply voltage and features an **analog interface**. The signal is internally processed and **temperature compensated**.

The outstanding performance of these sensors is based on Sensirion's **patented CMOSens® sensor technology**, which combines a thermal sensor element, signal processing and digital calibration on a single microchip. This assures **very fast signal processing**. An additional temperature sensor provides a highly **accurate temperature measurement in the gas stream**.

The well-proven CMOS technology is perfectly suited for high-quality mass production and is the ideal choice for demanding and cost-sensitive OEM applications.

Applications

- Ventilation
- Anesthesia
- Inspiratory flow sensing
- Gas mixing

OEM options

A variety of custom options can be implemented for highvolume OEM applications (digital interface, custom flow rates, calibration for other gases, different body form factor etc.). Contact us for more information.

Sensor electronics

The SFM3100 flow meter features a fifth-generation silicon sensor chip and provides an analog output signal. An EEPROM is included to store sensor specific information such as the serial number or customer specific data. See chapter 2.3 for more detail.

Direct gas temperature measurement

The SFM3100 is equipped with an NTC element centrally positioned in the main flow channel for direct temperature measurement of the gas.



1. Sensor Performance

1.1 Physical specifications ¹

Parameter	Condition	Vá	Value			
Flow Ranges	Air	-24	slm ²			
		Typ. ³	Max. ⁴			
Accuracy Flow 5,6,7	span (< 60 slm)	2.5	3	% m.v. ⁸		
	span (< 150 slm)	3.5	4	% m.v. ⁸		
	span (< 240 slm)	4.5	6	% m.v. ⁸		
	offset	0.02	0.03	slm ²		
Peopletion Flow (analog signal)	span	0.05	0.07	% m.v.		
Resolution Flow (analog signal)	offset	0.04	0.045	slm ²		
	span (<60 slm)	0.5	1	% m.v. ⁸		
	span (<150 slm)	1	2	% m.v. ⁸		
Noise Level 6,9	span (<240 slm)	2	3	% m.v. ⁸		
	offset	0.045	0.05	slm		
Accuracy Shift Due to Temperature	Span	0.4	0.5	% m.v. /10°C		
Variation ¹⁰	offset	0.015	0.02	slm/10°C		
Pressure Drop	@60slm	< 300 / < 1.2	< 350 / < 1.4			
	@200slm	< 1600 / < 6.4	< 1800 / < 7.2	Pa / inH₂O		
	0 °C - 80 °C		2.5	C°		
Accuracy Temperature	-20 °C - 0 °C		4	°C		
Resolution Temperature	analog signal	< 0.5°C ≈ 5mV/°C				

- ⁵ Including offset, non-linearity, hysteresis
- ⁶ Total accuracy/noise level/resolution is a sum of offset and span accuracy/noise level/resolution
- ⁷ The accuracy is valid for T(gas)=T(Chip)

¹Unless otherwise noted, all sensor specifications are valid at 25°C with Vdd = 5V and absolute pressure = 966 mbar and horizontal flow direction

² slm: mass flow measured in liters per minute at standard conditions (T = 20 °C, p = 1013.25 mbar)

³ for "Typ" a CpK of 0.67 is targeted (95% of sensors within the Typ limit)

⁴ for "Max" no sensor measured outside of this limits will be shipped and a CpK of 1.33 is targeted

⁸ %m.v. = % measured value = % of reading

⁹ noise level defined as standard deviation of individual sensor readings, measured at full sampling rate (typ: average of noise level ; max: at least 99.99% of sensors have a noise level below indicated value)

¹⁰ these effects need to be added to the initial values if applicable



1.2 Ambient conditions

Parameter Condition		Value	Unit
Calibrated Temperature Range	dry air	+5 +60	°C
Operating Temperature Range 10-95% rel. hum. (non cond.)		+5 +60	°C
Storage Temperature	0-100% rel. hum. (non cond.)	-25 +65	°C
Operating Pressure Range	absolute	0.6 – 1.07	bar
Burst Overpressure	gauge	>1000	mbar


1.3 Media compatibility

Parameter	Value
Calibration	Air
Media Compatibility	Air (non-condensing), N ₂ , O ₂ , other non-aggressive gases
Wetted Materials	PPE+PS blend, Si, Si $_3N_4$, SiO _x , Gold, Epoxy, Polyimide, silicone, stainless steel, FKM
RoHS, REACH	RoHS and REACH compliant



2. Electrical Specifications

2.1 Block diagram

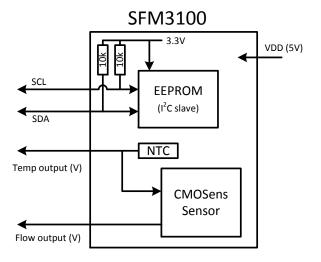

2.2 Electrical characteristics

Electrical properties	Condition		Value	Unit	
		Min.	Тур.	Max.	
Supply voltage (VDD)		4.75	5	5.25	V
Start-up Time ¹¹				200	ms
Power Consumption				150	mW
Flow output signal	range	0.095		2.45	V
Temperature output signal	range	0.2		2.4	V
EEPROM I ² C communication level	High	2.5	3.3	4.3	V
	Low	GND		0.8	V
Output capacitive load Cload				100	nF
Output resistive load R _{load}		10			kΩ

¹¹ After 4.75V is reached



2.3 EEPROM

The SFM3100 has an on-board EEPROM, which can be accessed over the I^2C interface.

The EEPROM is read-only and stores sensor specific information such as the serial number or customer specific data. For OEM options please contact Sensirion.


Address	Content
0xA002-0xA007	Product number (a-bbbbbb-cc)
	a: 0xA002 bbbbbb: 0xA003-0xA006 cc: 0xA007
0xA00A-0xA011	Serial number
	MSB: 0xA00A

The chip enable (E2) of the EEPROM is connected to ground so that the device address of the EEPROM starts with 10100b. The EEPROM is the M24C08-WDW6TP from ST Microelectronics. For more details please refer to the datasheet on the ST.com website.

2.4 Pin layout and mechanical concept of the electrical connection

Pin	Description	Pin	Description
1	Flow signal	2	Connect to GND
3	Temperature signal	4	GND
5	EEPROM SCL	6	EEPROM SDA
7	Not connected	8	VDD (5V)
9	GND	10	Not connected

The SFM3100 is designed for connector attachment and uses Molex Milli-Grid 87831-1041 (THT). Please refer to the datasheet of this part for mating connectors.

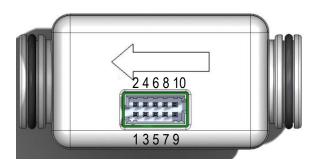

3. Conversion to Physical Values


3.1 Air flow calculation

In order to obtain the flow in slm from the voltage flow signal, the measured value (U[V]) needs to be converted using the following formula:


$$Q[slm] = \left(U[V] - \frac{Offset}{B}\right) \cdot \frac{A}{B} \cdot abs\left(U[V] - \frac{Offset}{B}\right)$$
$$= \left(U[V] - 0.67V\right) \cdot 80\frac{slm}{V^2} \cdot abs(U[V] - 0.67V)$$

The factors A, B and Offset are stored on the EEPROM. Please refer to Appendix B, formula (1).

For a flow smaller than -24.2 slm the sensor output is fixed at 0.095 V. Similarly, for flows above 240 slm the output voltage is set to 2.45 V.

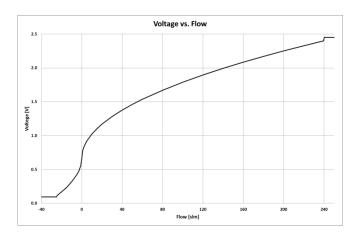

In order to obtain the measured temperature in °C, use the table attached in Appendix A or the following graph:

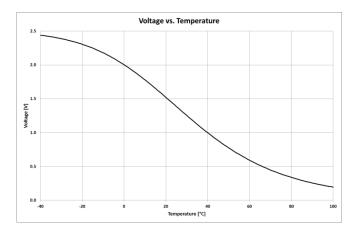

Please note that the first measurement performed directly after chip initialization is not valid.


3.2 O2 flow calculation


To measure oxygen flow a conversion function needs to be applied. Please refer to Appendix B.




4. Mechanical Specifications

All dimensions are in millimetres (mm).

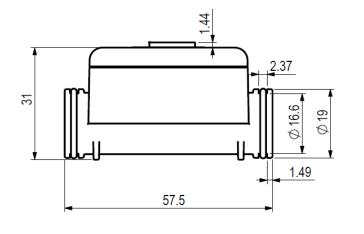

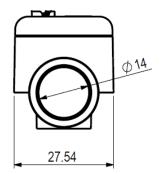

O-Ring NBR, ID 16.1x 1.6 mm 70 ±5 Shore A

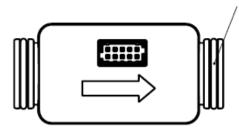




5. Instructions for Use

5.1 Calibration orientation

The sensors are calibrated horizontally as depicted in the following graph:

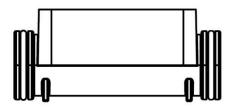

5.2 Inlet flow conditions

In order to provide good flow conditions, the inner diameter of the connecting tube has to be the same as the inner diameter of the SFM3100 main flow channel (16.6mm). The inlet tube is recommended to be straight and at least 10 cm in length. The SFM3100 is equipped with a mesh on the inlet of the flow channel to improve flow conditions.

Please refer to the application note "Inlet conditions for the SFM3000 Mass flow meters" for more information.

5.3 Sensor handling

The packaging method of the CMOSens chip together with the inert housing and the sealing materials ensure a tight and highly resistant sealing of the device. Please be aware that aggressive and corrosive gases can influence the sensor element and adversely affect the sealing and the plastic housing.

The use of explosive or toxic gases requires exceptional care and precautions as even the slightest leakage may result in dangerous situations.

For the above reasons, Sensirion only guarantees the safe use of the CMOSens® Mass Flow Meter for inert, inexplosive and non-toxic gases.

The SFM3100 sensor is designed to be robust and shock resistant. Nevertheless, the accuracy of the high-precision SFM3100 may be degraded when handled too rough. Sensirion does not guarantee proper operation in case of improper handling. **Note:** Please avoid applying mechanical stress to the sensor or the mesh.



The sensor is shipped in an antistatic tray to prevent damage by electrostatic discharge. In order to avoid damaging the sensor, ground yourself with a grounding strap or by touching a grounded object before touching the sensor. Furthermore it is advised to store the parts in an antistatic package when not in use.

5.4 ESD

The electronics of the SFM3100 sensor consist of an automotive qualified chip. It complies with the following ESD norms:

- AEC Q 100 002 (4kV HBM)
- AEC Q 100 003 (200V MM)

For ESD tests of the external EEPROM, please consult the datasheet of this component (see chapter 2.3).

Although the sensor complies with these norms, it does not mean the sensor is immune against ESD.

6. Ordering Information

Use the part names and product numbers shown in the table below when ordering SFM3100 sensors. For the latest product information and local distributors, visit <u>www.sensirion.com</u>.

Part name	Product Number 1-10xxxx-xx
SFM3100-240-VC	1-101290-02

Packaging units: 30 items/tray.

Every sensor is traceable by a unique Serial Number.

7. Revision history

Date	Version	Author	Changes		
Sept 2015	0.9	VRR	Preliminary Datasheet		
Dec 2015	0.91	VRR	Corrections:		
			 remove polyurethane from wetted materials list 		
			- adapt Appendix B		
Jan 2016	0.92	VRR	Reformulate formulas Appendix B, Chap. 3.1, add explanatory information		
May 2016	0.95	VRR	 improved accuracy specs, operating conditions 		
-			 revised EEPROM address of product number 		
			 change of product number 		
Dec 2016	1.0	SAW	Removal of preliminary, adjusted title, new part name SFM3100-240-VC		
Jun 2017	1.01	ALAN	Added Silicone to wetted materials		





Appendix A: Temperature Curve of additional temperature sensor

1.664		<u> </u>	Jerutui			aantio		iperat	ure sei	1301				
°C	kOhm	Volt	°C	kOhm	Volt	°C	kOhm	Volt	°C	kOhm	Volt	°C	kOhm	Volt
-40	1611	2.439	5	126.1	1.896	50	18.06	0.775	95	3.928	0.223	140	1.153	0.070
-39	1510	2.435	6	120.1	1.873	51	17.39	0.755	96	3.812	0.217	141	1.125	0.068
-38	1415	2.431	7	114.4	1.850	52	16.74	0.735	97	3.699	0.211	142	1.097	0.066
-37	1327	2.426	8	109.0	1.826	53	16.13	0.716	98	3.590	0.205	143	1.071	0.065
-36	1245	2.422	9	103.9	1.803	54	15.54	0.697	99	3.485	0.199	144	1.046	0.063
-35	1169	2.417	10	99.07	1.778	55	14.97	0.678	100	3.384	0.194	145	1.021	0.062
-34	1098	2.412	11	94.48	1.754	56	14.43	0.660	101	3.285	0.189	146	0.9965	0.060
-33	1031	2.406	12	90.12	1.729	57	13.91	0.643	102	3.191	0.184	147	0.9730	0.059
-32	969.3	2.400	13	85.99	1.704	58	13.41	0.625	103	3.099	0.179	148	0.9501	0.058
-31	911.3	2.394	14	82.07	1.678	59	12.93	0.608	104	3.010	0.174	149	0.9279	0.056
-30	857.2	2.388	15	78.35	1.652	60	12.47	0.592	105	2.925	0.170	150	0.9062	0.055
-29	806.6	2.381	16	74.82	1.626	61	12.03	0.576	106	2.842	0.165	151	0.8852	0.054
-28	759.3	2.374	17	71.47	1.600	62	11.61	0.560	107	2.762	0.161	152	0.8648	0.053
-27	715.0	2.367	18	68.29	1.574	63	11.20	0.545	108	2.684	0.156	153	0.8449	0.051
-26	673.6	2.359	19	65.26	1.547	64	10.81	0.530	109	2.609	0.152	154	0.8255	0.050
-25	634.8	2.351	20	62.39	1.520	65	10.44	0.515	110	2.536	0.148	155	0.8067	0.049
-24	598.5	2.343	21	59.65	1.493	66	10.08	0.501	111	2.466	0.144	156	0.7884	0.048
-23	564.5	2.334	22	57.06	1.467	67	9.735	0.487	112	2.398	0.141	157	0.7705	0.047
-22	532.6	2.325	23	54.59	1.440	68	9.404	0.474	113	2.332	0.137	158	0.7532	0.046
-21	502.7	2.315	24	52.23	1.413	69	9.085	0.461	114	2.269	0.134	159	0.7363	0.045
-20	474.6	2.305	25	50.00	1.386	70	8.778	0.448	115	2.207	0.130	160	0.7198	0.044
-19	448.3	2.294	26	47.87	1.359	71	8.484	0.436	116	2.147	0.127	161	0.7038	0.043
-18	423.6	2.283	27	45.84	1.332	72	8.201	0.424	117	2.090	0.124	162	0.6882	0.042
-17	400.4	2.272	28	43.91	1.305	73	7.928	0.412	118	2.034	0.120	163	0.6731	0.041
-16	378.6	2.260	29	42.07	1.278	74	7.666	0.400	119	1.979	0.117	164	0.6583	0.040
-15	358.1	2.248	30	40.32	1.252	75	7.414	0.389	120	1.927	0.114	165	0.6439	0.039
-14	338.8	2.235	31	38.65	1.225	76	7.171	0.378	121	1.876	0.111	166	0.6298	0.039
-13	320.7	2.222	32	37.06	1.199	77	6.938	0.368	122	1.826	0.109	167	0.6162	0.038
-12	303.7	2.208	33	35.54	1.173	78	6.713	0.358	123	1.778	0.106	168	0.6029	0.037
-11	287.6	2.193	34	34.09	1.147	79	6.496		124	1.732	0.103		0.5899	0.036
-10	272.5	2.179	35	32.71	1.122	80	6.288	0.338	125	1.687	0.101	170	0.5772	0.035
-9	258.3	2.163	36	31.39	1.096	81	6.087	0.329	126	1.643	0.098	171	0.5649	0.035
-8	244.9	2.147	37	30.13	1.071	82	5.894	0.320	127	1.601	0.096	172	0.5529	0.034
-7	232.3	2.131	38	28.93	1.046	83	5.707	0.311	128	1.560	0.093	173	0.5412	0.033
-6	220.4	2.114	39	27.78	1.022	84	5.528	0.302	129	1.520	0.091	174	0.5298	0.033
-5	209.2	2.097	40	26.69	0.998	85	5.354	0.294	130	1.482	0.089	175	0.5187	0.032
-4	198.6	2.079	41	25.64	0.974	86		0.286	131	1.444	0.087	176	0.5078	0.031
-3	188.6	2.061	42	24.64	0.950	87	5.027	0.278	132	1.408	0.085	177		0.031
-2	179.1	2.042	43	23.68	0.927	88		0.270	133	1.372	0.083	178	0.4869	0.030
-1	170.2	2.022	44	22.77	0.904	89	4.722	0.263	134	1.338	0.081		0.4769	0.029
0	161.8	2.002	45	21.90	0.882	90	4.578	0.256	135	1.305	0.079	180	0.4670	0.029
1	153.9	1.982	46	21.06	0.859	91	4.438	0.249	136	1.273	0.077			
2	146.3	1.961	47	20.26	0.838	92	4.304	0.242	137	1.241	0.075			
3	139.2	1.940	48	19.49	0.816	93			138	1.211	0.073			
4	132.5	1.918	49	18.76	0.795	94	4.049	0.229	139	1.181	0.071			



Appendix B: Oxygen flow conversion

The O2 flow can be calculated using the following two formulas.

(1)
$$Q[slm] = \left(U[V] - \frac{Offset}{B}\right) \cdot \frac{A}{B} \cdot abs\left(U[V] - \frac{Offset}{B}\right)$$

(2)
$$Q_{02}[slm] = \frac{Q[slm]}{\alpha} \left\{ 1 - \frac{F_{02,i}}{\beta} \left(\frac{C_{02} - 21\%}{79\%} \right) \right\}$$

U[V] = flow signal from sensor (voltage)

 C_{02} = O₂ concentration (e.g. 100%, if pure O₂ is used)

Q[slm] = calculation of flow without O₂ correction

refer to chapter 3.1: flow calculation

 $Q_{02}[slm]$ = flow at a specified O₂ concentration

All other constants for formula (1) and (2) and the matrix for the O_2 correction values ($Q_{O2} \times F_{O2}$) can be found on the EEPROM. The values are stored as 16bit Int.

Use value $F_{O2,i}$ where $Q_{O2,i}$ is closest to Q[slm] (from (1)).

EEPROM addresses

Value	HEX Address	Description	
A	0xA1AA	Boguired for formula (1)	
A	0xA1AB	Required for formula (1)	
В	0xA1A8	Deguired for formula (1)	
D	0xA1A9	Required for formula (1)	
Offset	0xA1AC	Deguired for formula (1)	
Olisei	0xA1AD	Required for formula (1)	
<i></i>	0xA100	Deguired for formula (2)	
α	0xA101	Required for formula (2)	
0	0xA102	Boguired for formula (2)	
β	0xA103	Required for formula (2)	

Q ₀₂ x F ₀₂						
size of matrix	0xA106	Size of the Qo2 x Fo2 matrix				
SIZE OF MALIIX	0xA107					
0	0xA108	Lowest flow, e.g40ml/s (value as 16bit Int)				
Q _{02,1}	0xA109	required for formula (2)				
r	0xA10A	Correction factor corresponding to Q _{02,1} , e.g. 0.02				
F _{02,1}	0xA10B	(value as 16bit Int) required for formula (2)				
	0xA1	Bequired for formula (2)				
	0xA1	Required for formula (2)				
	0xA1	Dequired for formula (2)				
	0xA1	Required for formula (2)				
0	0xA1A4	Dequired for formula (2)				
Q _{02,40}	0xA1A5	Required for formula (2)				
r	0xA1A6	Deswired for formula (2)				
F _{02,40}	0xA1A7	Required for formula (2)				



Important Notices

Warning, personal injury

Do not use this product as safety or emergency stop devices or in any other application where failure of the product could result in personal injury (including death). Do not use this product for applications other than its intended and authorized use. Before installing, handling, using or servicing this product, please consult the datasheet and application notes. Failure to comply with these instructions could result in death or serious injury.

If the Buyer shall purchase or use SENSIRION products for any unintended or unauthorized application, Buyer shall defend, indemnify and hold harmless SENSIRION and its officers, employees, subsidiaries, affiliates and distributors against all claims, costs, damages and expenses, and reasonable attorney fees arising out of, directly or indirectly, any claim of personal injury or death associated with such unintended or unauthorized use, even if SENSIRION shall be allegedly negligent with respect to the design or the manufacture of the product.

ESD Precautions

The inherent design of this component causes it to be sensitive to electrostatic discharge (ESD). To prevent ESD-induced damage and/or degradation, take customary and statutory ESD precautions when handling this product.

See application note "Handling Instructions" for more information.

Warranty

SENSIRION warrants solely to the original purchaser of this product for a period of 12 months (one year) from the date of delivery that this product shall be of the quality, material and workmanship defined in SENSIRION's published specifications of the product. Within such period, if proven to be defective, SENSIRION shall repair and/or replace this product, in SENSIRION's discretion, free of charge to the Buyer, provided that:

notice in writing describing the defects shall be given to SENSIRION within fourteen (14) days after their appearance;

Headquarters and Subsidiaries

SENSIRION AG Laubisruetistr. 50 CH-8712 Staefa ZH Switzerland

phone: +41 44 306 40 00 fax: +41 44 306 40 30 info@sensirion.com www.sensirion.com

Sensirion AG (Germany) phone: +41 44 927 11 66 info@sensirion.com www.sensirion.com Sensirion Inc., USA phone: +1 805 409 4900 info_us@sensirion.com www.sensirion.com

Sensirion Japan Co. Ltd. phone: +81 3 3444 4940 info@sensirion.co.jp www.sensirion.co.jp such defects shall be found, to SENSIRION's reasonable satisfaction, to have arisen from SENSIRION's faulty design, material, or workmanship;

- the defective product shall be returned to SENSIRION's factory at the Buyer's expense; and
- the warranty period for any repaired or replaced product shall be limited to the unexpired portion of the original period.

This warranty does not apply to any equipment which has not been installed and used within the specifications recommended by SENSIRION for the intended and proper use of the equipment. EXCEPT FOR THE WARRANTIES EXPRESSLY SET FORTH HEREIN, SENSIRION MAKES NO WARRANTIES, EITHER EXPRESS OR IMPLIED, WITH RESPECT TO THE PRODUCT. ANY AND ALL WARRANTIES, INCLUDING WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE EXPRESSLY EXCLUDED AND DECLINED.

SENSIRION is only liable for defects of this product arising under the conditions of operation provided for in the datasheet and proper use of the goods. SENSIRION explicitly disclaims all warranties, express or implied, for any period during which the goods are operated or stored not in accordance with the technical specifications.

SENSIRION does not assume any liability arising out of any application or use of any product or circuit and specifically disclaims any and all liability, including without limitation consequential or incidental damages. All operating parameters, including without limitation recommended parameters, must be validated for each customer's applications by customer's technical experts. Recommended parameters can and do vary in different applications.

SENSIRION reserves the right, without further notice, (i) to change the product specifications and/or the information in this document and (ii) to improve reliability, functions and design of this product.

Copyright © 2001-2015, SENSIRION. CMOSens[®] is a trademark of Sensirion All rights reserved

> Sensirion Korea Co. Ltd. phone: +82 31 337 7700~3 info-kr@sensirion.com www.sensirion.co.kr

Sensirion China Co. Ltd. phone: +86 755 8252 1501 info@sensirion.com.cn www.sensirion.com.cn

To find your local representative, please visit www.sensirion.com/contact